

7-LIGHT EMS PHOTON REJUVENATION VIBRATING BEAUTY DEVICE



NEGATIVE IONS
FAR INFRA RED (F.I.R.)
GERMANIUM
NANO SILVER

The first and only Beauty Device
With 7 types of therapies and
4 additional elements



- 7 therapy selection !**
- >>FIRMING
 - >>LIFTING
 - >>SMOOTH
 - >>DELICATE
 - >>RADIANT COMPLEXION



7-LIGHT THERAPY



Red Light:

- Stimulates collagen production, improves blood circulation, and promotes healing.
- Reduces wrinkles and fine lines, improves skin texture and tone, promotes overall rejuvenation and anti-aging effects.



Blue Light:

- Kills acne-causing bacteria, reduces inflammation, and regulates oil production.
- Controls acne breakouts, treats acne scars and blemishes, and helps minimize pores.



Green Light:

- Calms and soothes the skin, reduces redness and irritation.
- Balances skin tone, reduces hyperpigmentation and sunspots, provides a calming effect on sensitive skin.



Yellow Light:

- Boosts lymphatic circulation, improves cellular oxygen exchange.
- Promotes detoxification of the skin, enhances blood flow, reduces redness and inflammation.



7-LIGHT THERAPY



Orange Light:

- Increases energy in skin cells, revitalizes dull complexion.
- Restores radiance and brightens the skin, improves overall skin tone.



Purple Light:

- Combines the benefits of red and blue light.
- Treats acne, promotes skin healing, improves overall complexion.



Cerulean Blue Light:

- Calms the skin, reduces inflammation and sensitivity.
- Soothes irritated skin, reduces redness, and helps alleviate skin conditions such as eczema and rosacea.

NOTE : *These functions and benefits may vary depending on the specific device being used and individual skin conditions.




7-LIGHT EMS PHOTON REJUVENATION VIBRATING BEAUTY DEVICE

The EMS 7-Light Photon Skin Rejuvenation Device offers several benefits, including:

- **Improved skin health:** The device's seven different light colors target various skin concerns, such as acne, wrinkles, pigmentation, and dullness. Scientifically supported photon therapy helps to stimulate collagen production, improve blood circulation, and rejuvenate the skin, resulting in a more radiant complexion.
- **Enhanced skincare product absorption:** The device's EMS (Electrical Muscle Stimulation) technology helps to enhance the absorption of skincare products. When used in conjunction with serums or creams, the EMS function aids in delivering active ingredients deep into the skin, maximizing their effectiveness.
- **Minimized appearance of wrinkles and fine lines:** The stimulation provided by EMS and specific light wavelengths helps to boost collagen and elastin production, which are essential for maintaining skin elasticity and firmness. This can lead to a reduction in the appearance of wrinkles and fine lines over time.
- **Reduced acne and inflammation:** The blue light emitted by the device has antibacterial properties that can target acne-causing bacteria and help reduce inflammation. It can also regulate sebum production, preventing clogged pores and breakouts.
- **Improved skin tone and texture:** The device's various light colors, such as red, green, and yellow, can help even out skin tone, reduce hyperpigmentation, and enhance overall skin texture. This can result in a smoother and more uniform complexion.
- **Non-invasive and convenient:** The EMS 7-Light Photon Skin Rejuvenation Device offers a non-invasive and convenient at-home skincare solution. It allows users to customize their treatment according to their specific skincare needs and can be easily incorporated into their skincare routine.



ADDITIONAL 4 ELEMENTS

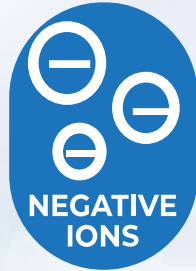


Negative Ions:

- Negative ions are molecules that have gained an extra electron, and they are found abundantly in nature, especially near waterfalls, forests, and beaches. When applied to the skin, negative ions help improve skin hydration and moisture retention.
- Increased skin hydration, improved skin elasticity, reduction in fine lines and wrinkles, enhanced absorption of skincare products.




Far Infrared Red:

- Far infrared red waves penetrate deeply into the skin, promoting blood circulation and stimulating cellular metabolism.
- Increased collagen production, improved skin texture and tone, reduced appearance of scars and pigmentation, enhanced skin rejuvenation and healing.

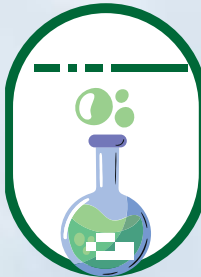

Germanium:

- Germanium is a natural element that emits far infrared rays and negative ions. When applied to the skin, Germanium helps increase blood flow and oxygen supply to the cells.
- Improved skin vitality, boosted collagen production, increased skin elasticity, reduced signs of aging, enhanced detoxification and cleansing of the skin.

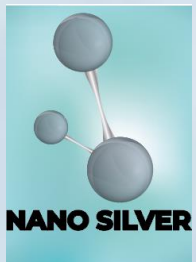

Nano Silver:

- Nano Silver particles possess antimicrobial properties, effectively fighting against harmful bacteria and reducing inflammation on the skin.
- Reduced acne and breakouts, improved wound healing, minimized skin irritations and redness, enhanced overall skin health and cleanliness.





WITH EMS TECHNOLOGY

EMS stands for Electrical Muscle Stimulation. In a light therapy device, EMS refers to the use of low-level electrical currents to stimulate the facial muscles.

1. Lifts and Tightens: EMS helps to lift and tighten sagging facial muscles, giving a more youthful and lifted appearance to the skin.

2. Defines Facial Contours: By toning the muscles, EMS can help define the jawline and cheekbones, resulting in a more sculpted and defined facial appearance.

3. Reduces Wrinkles and Fine Lines: As the muscles are stimulated and strengthened, the appearance of wrinkles and fine lines can be reduced.

4. Improves Circulation: EMS stimulates blood circulation in the skin, bringing more oxygen and nutrients to the cells. This can result in improved skin health and a natural glow.

5. Enhances Product Absorption: The muscle contractions caused by EMS can help in the better absorption of skincare products into the deeper layers of the skin. This maximizes the effectiveness of the products used in combination with the device.

6. Relaxation and Stress Relief: EMS can provide a soothing and relaxing effect on the facial muscles, helping to relieve tension and stress in the face.



WITH VIBRATION EFFECT



The EMS 7-Light Photon Skin Rejuvenation Vibrating Device utilizes a vibrating effect as part of its functionality. This vibrating effect helps to stimulate and massage the skin, promoting blood circulation and lymphatic drainage.

The vibrations also facilitate deeper penetration of skincare products into the skin, enhancing their absorption and maximizing their benefits. Furthermore, the gentle vibrations can provide a soothing and relaxing sensation during use, adding to the overall user experience.



HOW TO USE



1. Cleanse your face: Start by cleansing your face thoroughly to remove any dirt, makeup, or skincare products.
2. Prep your skin: Apply any serums or toners that you usually use before starting your skincare routine. Make sure your skin is clean and dry.
3. Choose the light setting: Light therapy devices usually come with different light settings, such as red, blue, etc. Select the light setting that best suits your skincare needs.
4. Turn on the device: Switch on the light therapy device and ensure it is fully charged.
5. Start the treatment: Gently place the device on your skin and move it in circular motions or go over specific areas.
6. Be consistent: For effective results, it's essential to use the light therapy device consistently. Follow the specific instructions provided with your device regarding frequency and duration of treatments.
7. Post-treatment skincare: After your light therapy session, you can apply your usual moisturizer or any post-treatment skincare products you typically use. This will help lock in the benefits of the light therapy treatment.




HOW TO USE

8. Clean and store the device: Clean the light therapy device as per the manufacturer's instructions after each use.

Make sure to store it in a safe place, away from excessive heat or moisture.

9. Duration of use: Light therapy sessions can range from 10 minutes to around 20 minutes

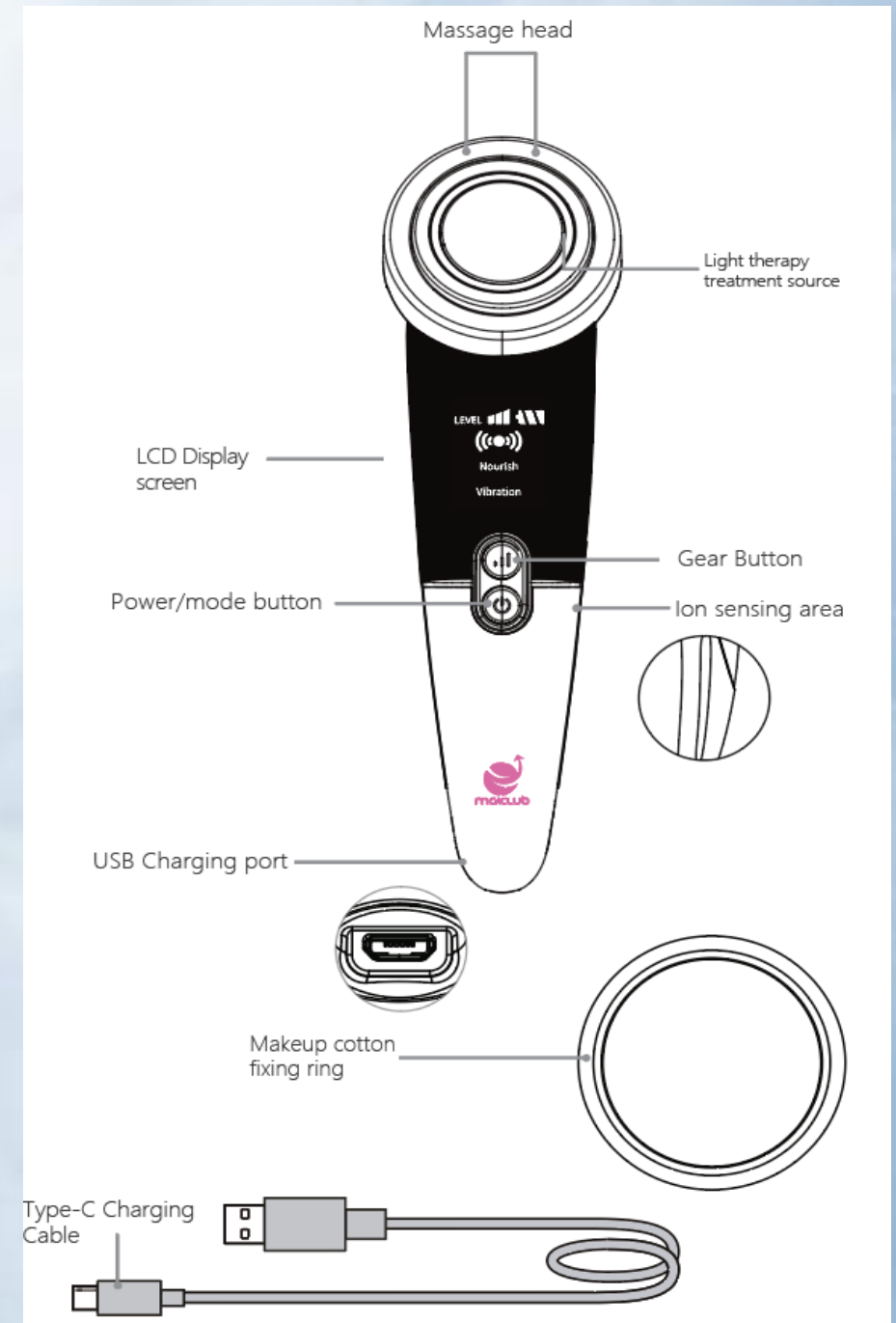



FREQUENTLY ASKED QUESTIONS (FAQ)

Question	Solution
Why is the massage head not vibrating?	Hold the beauty device with your fingers on both sides. This allows the massage head to come into contact with the skin and form a circuit with the body.
Why doesn't it vibrate when used by someone else?	If you want to use it on someone else's face, the device needs to have physical contact with the user to form a circuit. Otherwise, the device will not function properly.
How to clean it after use?	After use, dampen a clean piece of cotton and gently wipe the head of the device.
How long does it take to charge?	Charging time is approximately 3 hours. The device cannot be turned on while charging, so please use it after charging.

7-EMS BEAUTY DEVICE CE REPORT



T: 0086-755-2523 4088
E-mail: agc@agccert.com
Web: Http://www.agccert.com
1-2/F, Building 19, Junfeng Industrial Park
Chongqing Road, Heping Community, Fuhai Street
Bao'an District, Shenzhen, Guangdong, China

Attestation of Conformity

Registration No. AGC15789240602E2

Applicant MEGA GLORYOUNG INTERNATIONAL CORP
Rm 315, Filamco Building, Soler Street Corner Masangkay Street,
Binondo, Manila, Philippines.

Product Designation EMS 7-LIGHT PHOTON SKIN REJUVENATION VIBRATING BEAUTY DEVCE

Brand Name MGICLUB

Model / Series Models MGI-080

Requirement	Applied Standards	Document Evidence	Result
EMC Directive	EN IEC 55014-1:2021 EN IEC 55014-2:2021 EN IEC 61000-3-2:2019/A1:2021 EN 61000-3-3:2013/A2:2021	Test Report: AGC15789240602EE01	Conform



Signed by General Manager(King Zhang)
Issue Date: July 8, 2024



This Attestation of Conformity is recognized by Attestation of Global Compliance (Shenzhen) Co., Ltd. and made in accordance with the EMC Directive 2014/30/EU. The attestation doesn't imply assessment of the production. The Applicant of this attestation is authorized to use this attestation in connection with EC declaration of conformity to the Directive. The attestation is only applicable to the equipments described above. This attestation shall not be re-produced except in full without the written approval of Attestation of Global Compliance (Shenzhen) Co., Ltd.

Note: This attestation is part of the full test report(s) and should be used in conjunction with it.



Recognized by China National Accreditation Service for Conformity Assessment

T: 0086-755-2523 4088
E-mail: agc@agccert.com
Web: www.agccert.com
6/F., Building 2, Sanwei Chaxi Industrial Park, Sanwei Community, Hangcheng Street, Bao'an District, Shenzhen, Guangdong, China

Attestation of Conformity

Registration No.: AGC15789240602-C001

Holder: Name: MEGA GLORYOUNG INTERNATIONAL CORP
Address: Rm 315, Filamco Building, Soler Street Corner Masangkay Street,
Binondo, Manila, Philippines.

Product Description of the Attestation: Product Name: EMS 7-LIGHT PHOTON SKIN REJUVENATION VIBRATING BEAUTY DEVCE
Brand Name: MGICLUB
Model / Series models: MGI-080

This attestation of conformity is according to the provisions of EU ROHS Directive 2011/65/EU& (EU) 2015/863, it is confirmed the product has been tested and found in conformity with below standards.

Requirement	Applied Standards	Document Evidence	Result
2011/65/EU & (EU) 2015/863	IEC 62321-3-1:2013 / IEC 62321-5:2013	Test Report: AGC15789240602-001	Conform
	IEC 62321-4:2013+A1:2017		
	IEC 62321-7-2:2017/ IEC 62321-6:2015		
	IEC 62321-7-1:2015/ IEC 62321-8:2017		



Signed by: Technical Director (Leon)
Issued Date: Jun. 27, 2024








The attestation is only applicable to the above product submitted by applicant. It doesn't imply assessment of the whole production. The applicant is authorized to use this attestation in connection with attestation of conformity to Directive 2011/65/EU & (EU) 2015/863. This attestation may not be partially reproduced without the approval of Attestation of Global Compliance (Shenzhen) Std & Tech Co., Ltd.

Note: This declaration is part of the full test report(s) and should be used in conjunction with it, and if the report(s) is invalid then the declaration is invalid.

7-EMS BEAUTY DEVICE CE REPORT

AGC®

Report No.: AGC15789240602-001

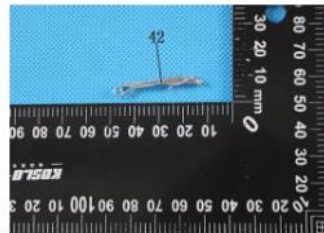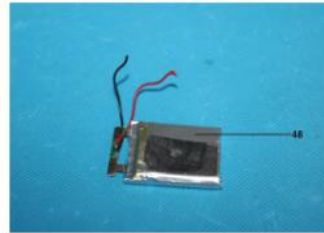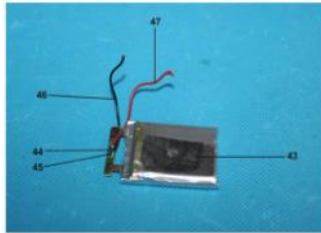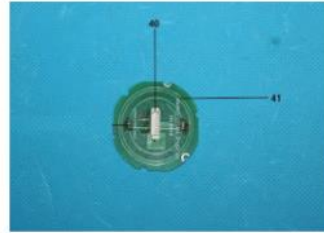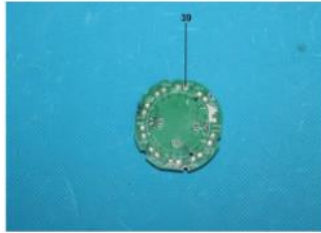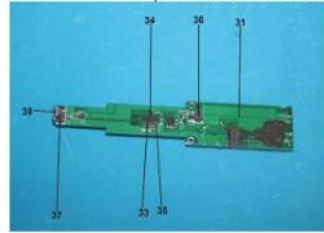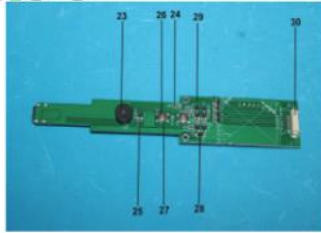

Any report having not been signed by authorized approver, or having been altered without authorization, or having not been stamped by the "Dedicated Testing/Inspection Stamp" is deemed to be invalid. Copying or excerpting portion of, or altering the content of the report is not permitted without the written authorization of AGC. The test results presented in the report apply only to the tested sample. Any objections to report issued by AGC should be submitted to AGC within 15days after the issuance of the test report. Further enquiry of validity or verification of the test report should be addressed to AGC by agc01@agccert.com.

Attestation of Global Compliance(Shenzhen)Co., Ltd

AGC®

Report No.: AGC15789240602-001



The photo of AGC15789240602-001 is for use only with the original report.

Test Point Description

Test point	Test module	Test parts	Test point description
EMS 7-LIGHT PHOTON SKIN REJUVENATION VIBRATING BEAUTY DEVCE Model : MGI-080			
1	---	---	Translucent plastic ring
2	---	---	Plastic ring with silver electroplating
3	---	---	Milk white plastic ring
4	---	---	Silver metallic sheet
5	---	---	Milk white sealing ring
6	---	---	Khaki cloth
7	---	---	Grey white silicone sheet
8	---	---	Brown transparent plastic shell
9	---	---	Silver coating
10	---	---	White plastic shell
11	---	---	Champagne coating
12	---	---	Black foam with glue
13	---	---	Metal weight block
14	---	---	Silver screw
15	---	---	Metallic shell
16	---	Motor	Red wire jacket
17	---	---	Blue wire jacket
18	---	Terminal wire	Red gray wire jacket
19	Circuit board	Digital tube	Black plastic sheet
20			White plastic shell
21			PCB
22			Milk white glue
23	---	Buzzer	Black plastic shell
24	---	---	Chip capacitor
25	---	---	Chip resistor
26	Key	---	Copper metal button
27			Metallic shell
28	---	---	Chip diode

Any report having not been signed by authorized approver, or having been altered without authorization, or having not been stamped by the "Dedicated Testing/Inspection Stamp" is deemed to be invalid. Copying or excerpting portion of, or altering the content of the report is not permitted without the written authorization of AGC. The test results presented in the report apply only to the tested sample. Any objections to report issued by AGC should be submitted to AGC within 15days after the issuance of the test report. Further enquiry of validity or verification of the test report should be addressed to AGC by agc01@agccert.com.

Attestation of Global Compliance(Shenzhen)Co., Ltd

7-EMS BEAUTY DEVICE CE REPORT




Report No.: AGC15789240602EE01
Page 14 of 50

5.5. Test Summary

Test Engineer	Jimu Lao	Temperature	22.5 °C
Test Date	2024-06-15	Air Pressure	985 Mbar
Worst Mode	Mode 2	Relative Humidity	68.8 %
Verdict	Pass		

Test graph and data for Conducted Emission at line L1

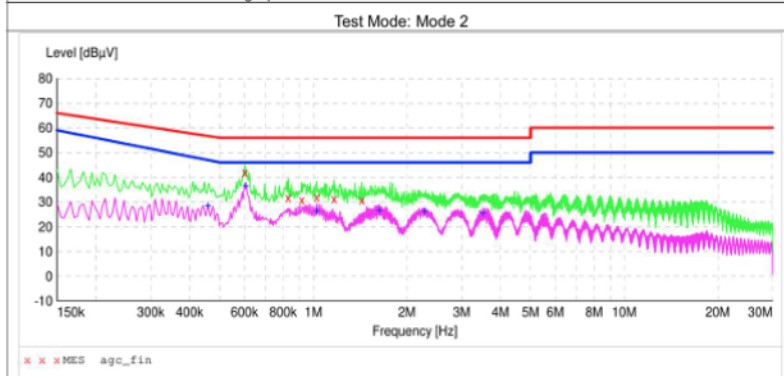

Frequency[MHz]	Level[dBμV]	Factor[dB]	Limit[dBμV]	Margin[dB]	Detector	Line
0.602000	41.9	6.2	56.0	14.1	QP	L1
0.830000	31.9	6.2	56.0	24.1	QP	L1
0.918000	30.9	6.2	56.0	25.1	QP	L1
1.026000	32.1	6.2	56.0	23.9	QP	L1
1.166000	31.3	6.2	56.0	24.7	QP	L1
1.434000	30.8	6.2	56.0	25.2	QP	L1
0.458000	29.1	6.1	46.9	17.8	AV	L1
0.606000	37.1	6.2	46.0	8.9	AV	L1
1.026000	27.0	6.2	46.0	19.0	AV	L1
1.630000	27.2	6.2	46.0	18.8	AV	L1
2.282000	26.8	6.3	46.0	19.2	AV	L1
3.518000	26.2	6.3	46.0	19.8	AV	L1


Any report having not been signed by authorized approver, or having been altered without authorization, or having not been stamped by the "Dedicated Testing/Inspection Stamp" is deemed to be invalid. Copying or excerpting portion of, or altering the content of the report is not permitted without the written authorization of AGC. The test results presented in the report apply only to the tested sample. Any objections to report issued by AGC should be submitted to AGC within 15days after the issuance of the test report. Further enquiry of validity or verification of the test report should be addressed to AGC by agc01@agccert.com.

Attestation of Global Compliance(Shenzhen)Co., Ltd
Attestation of Global Compliance(Shenzhen)Std & Tech Co., Ltd
Tel: +86-755 2523 4088 E-mail: agc@agccert.com Web: http://www.agccert.com/



Report No.: AGC15789240602EE01
Page 15 of 50

Test graph and data for Conducted Emission at line N

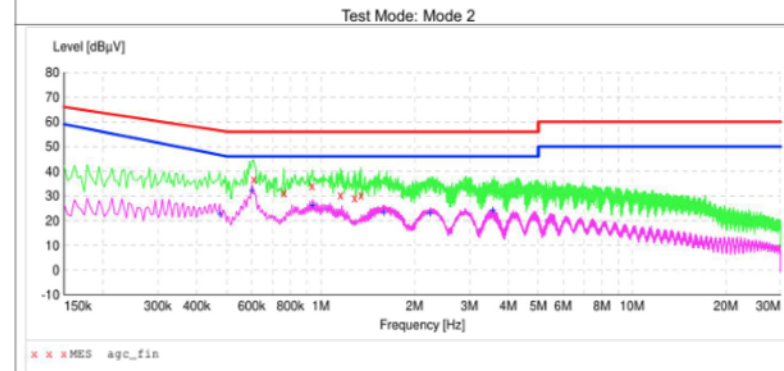

Frequency[MHz]	Level[dBμV]	Factor[dB]	Limit[dBμV]	Margin[dB]	Detector	Line
0.610000	36.8	6.2	56.0	19.2	QP	N
0.762000	31.3	6.2	56.0	24.7	QP	N
0.938000	34.2	6.2	56.0	21.8	QP	N
1.158000	30.4	6.2	56.0	25.6	QP	N
1.286000	29.2	6.2	56.0	26.8	QP	N
1.346000	30.5	6.2	56.0	25.5	QP	N
0.474000	23.2	6.1	46.6	23.4	AV	N
0.602000	32.8	6.2	46.0	13.2	AV	N
0.942000	26.8	6.2	46.0	19.2	AV	N
1.594000	23.7	6.2	46.0	22.3	AV	N
2.246000	23.4	6.3	46.0	22.6	AV	N
3.566000	24.8	6.3	46.0	21.2	AV	N

Any report having not been signed by authorized approver, or having been altered without authorization, or having not been stamped by the "Dedicated Testing/Inspection Stamp" is deemed to be invalid. Copying or excerpting portion of, or altering the content of the report is not permitted without the written authorization of AGC. The test results presented in the report apply only to the tested sample. Any objections to report issued by AGC should be submitted to AGC within 15days after the issuance of the test report. Further enquiry of validity or verification of the test report should be addressed to AGC by agc01@agccert.com.

Attestation of Global Compliance(Shenzhen)Co., Ltd
Attestation of Global Compliance(Shenzhen)Std & Tech Co., Ltd
Tel: +86-755 2523 4088 E-mail: agc@agccert.com Web: http://www.agccert.com/

7-EMS BEAUTY DEVICE CE REPORT




Report No.: AGC15789240602EE01
Page 44 of 50

15. Photographs of EUT



All view of EUT

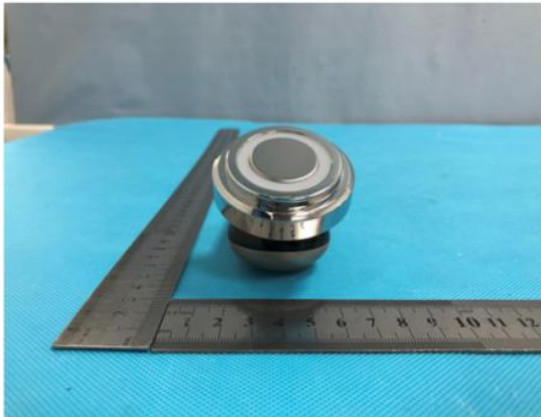

Top view of EUT

Any report having not been signed by authorized approver, or having been altered without authorization, or having not been stamped by the "Dedicated Testing/Inspection Stamp" is deemed to be invalid. Copying or excerpting portion of, or altering the content of the report is not permitted without the written authorization of AGC. The test results presented in the report apply only to the tested sample. Any objections to report issued by AGC should be submitted to AGC within 15days after the issuance of the test report. Further enquiry of validity or verification of the test report should be addressed to AGC by agc01@agccert.com.

Attestation of Global Compliance(Shenzhen)Co., Ltd
Attestation of Global Compliance(Shenzhen)Std & Tech Co., Ltd
Tel: +86-755 2523 4088 E-mail: agc@agccert.com Web: <http://www.agccert.com/>



Report No.: AGC15789240602EE01
Page 47 of 50



Right view of EUT



Open view of EUT


Any report having not been signed by authorized approver, or having been altered without authorization, or having not been stamped by the "Dedicated Testing/Inspection Stamp" is deemed to be invalid. Copying or excerpting portion of, or altering the content of the report is not permitted without the written authorization of AGC. The test results presented in the report apply only to the tested sample. Any objections to report issued by AGC should be submitted to AGC within 15days after the issuance of the test report. Further enquiry of validity or verification of the test report should be addressed to AGC by agc01@agccert.com.

Attestation of Global Compliance(Shenzhen)Co., Ltd
Attestation of Global Compliance(Shenzhen)Std & Tech Co., Ltd
Tel: +86-755 2523 4088 E-mail: agc@agccert.com Web: <http://www.agccert.com/>


7-EMS BEAUTY DEVICE CE REPORT



Report No.: AGC15789240602EE01
Page 49 of 50

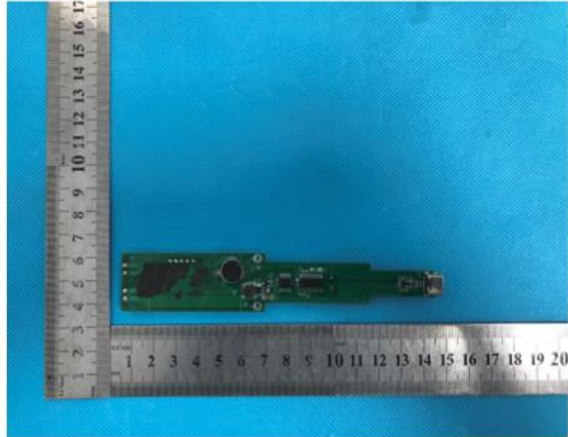

Internal view-3 of EUT

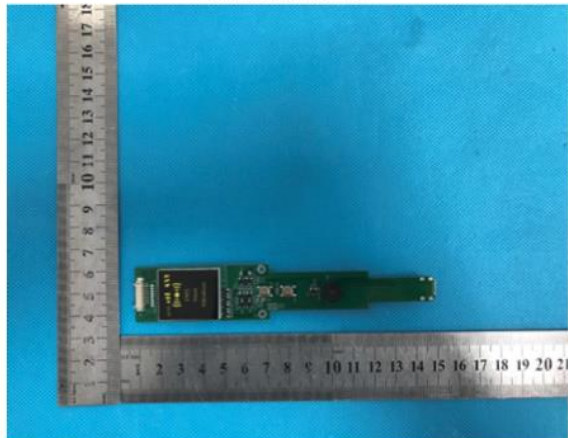

Internal view-4 of EUT

Any report having not been signed by authorized approver, or having been altered without authorization, or having not been stamped by the "Dedicated Testing/Inspection Stamp" is deemed to be invalid. Copying or excerpting portion of, or altering the content of the report is not permitted without the written authorization of AGC. The test results presented in the report apply only to the tested sample. Any objections to report issued by AGC should be submitted to AGC within 15days after the issuance of the test report. Further enquiry of validity or verification of the test report should be addressed to AGC by agc01@agccert.com.

Attestation of Global Compliance(Shenzhen)Co., Ltd
Attestation of Global Compliance(Shenzhen)Std & Tech Co., Ltd
Tel: +86-755 2523 4088 E-mail: agc@agccert.com Web: <http://www.agccert.com/>



Report No.: AGC15789240602EE01
Page 50 of 50



View of Battery

---End of Report---

Any report having not been signed by authorized approver, or having been altered without authorization, or having not been stamped by the "Dedicated Testing/Inspection Stamp" is deemed to be invalid. Copying or excerpting portion of, or altering the content of the report is not permitted without the written authorization of AGC. The test results presented in the report apply only to the tested sample. Any objections to report issued by AGC should be submitted to AGC within 15days after the issuance of the test report. Further enquiry of validity or verification of the test report should be addressed to AGC by agc01@agccert.com.

Attestation of Global Compliance(Shenzhen)Co., Ltd
Attestation of Global Compliance(Shenzhen)Std & Tech Co., Ltd
Tel: +86-755 2523 4088 E-mail: agc@agccert.com Web: <http://www.agccert.com/>



DISCLAIMER

It is important to note that EMS should be used in accordance with the device's instructions and should not be used excessively or without proper guidance.

It is recommended to consult with a skincare professional before starting any new skincare routine or using an EMS device.